

SOFTWARE GUIDE

INTRACRANIAL PRESSURE WAVEFORM GENERATOR

April 16, 2019

Senior Design Team 3
Wright State University
Department of Biomedical, Industrial, and Human Factors Engineering

Contents

Raspberry Pi Setup	4
Troubleshooting	6
Digital Analog Converting Software	7
Frequency Calculation	7
Computer Application	8
Waveform Codes	14

List of Figures

1	GUI	8
2	Normal Waveform	15
3	Abnormal Waveform	18
4	Leveled Waveform	21

Listings

1	Gui Code	9
2	Calibration	14
3	Normal Waveform	16
4	Abnormal Waveform	19
5	Leveled Waveform	22

RASPBERRY PI SETUP

1. Hardware

Raspberry Pi 3

Pi Cobbler

Adafruit Digital Analog Converter (DAC)

32 GB or more microSD Card

Keyboard

Computer Mouse

Wireless HDMI

2. Software Installation

NOOBS Install

We installed NOOBS as the operating system installer. A microSD card with at least 32 GB of storage is required. The set up was as follows:

1. Insert microSD into computer
2. Go to <https://www.raspberrypi.org/downloads/noobs/> and download the NOOBS zip file
3. Extract the zip file onto the microSD
4. Eject the microSD and insert into the Raspberry Pi 5. Power the Raspberry Pi on. The first screen that will appear select Raspbian and install
6. The Pi will now reboot and after it has power back on the Raspberry Pi is now ready for use

Set up: <https://www.youtube.com/watch?v=vvxCNQ5AYPg>

3. Login Credentials

The username is "pi" – all lowercase

The password is "ICP" – all capitalized letters

4. Package Installation

The process for installing any python packages should be done as follows from inside the terminal:

1. Run "sudo apt update" through the terminal before installing new software packages. This will ensure that all packages are the latest version and avoid any versioning issues with newer installed packages.
2. Type "sudo pip3 install xxxxxx" where xxxxxx is the name of the package with which you wish to install.
3. This will not include all necessary packages so if one of them doesn't work using "sudo pip3 install" then it is acceptable to use "sudo apt-get install python3-xxxxxx".

Viewing installed packages

To check if a package is installed type into the terminal "dpkg -l | grep <keyword>", where <keyword> is (part of) the name of the package you are interested in finding. For example, searching "dpkg -l | grep nvidia" will list all installed packages with "nvidia" in the name or description. This method will not work for meta-packages or repositories, though it will work for most cases. You can check where a package is installed by typing "sudo dpkg -S packagename".

The packages that are required for our system are as follows:

1. Adafruit MCP4275

Used for the DAC chip. Install was set up using:

<https://learn.adafruit.com/mcp4725-12-bit-dac-with-raspberry-pi>

2. Guizero

Used for the creation of the Graphical User Interface (GUI).

3. SciPy

contains modules for optimization, linear algebra, integration, interpolation, special functions, FFT, signal and image processing, ODE solvers and other tasks.

4. Numpy

Used for multi-dimensional arrays and matrices, along with a large collection of high-level mathematical functions.

5.I2C

A bus that allows a connection between the pi and wiring through the GPIO pins. Setup was

done through:

<https://learn.sparkfun.com/tutorials/raspberry-pi-spi-and-i2c-tutorial/all>

6. GPIO

These are the pins that the connects the code to the DAC chip. This connection is done through the Pi Cobbler. This software allows you to pick the pin the signal will come out of. Setup for this was done through:

<https://www.raspberrypi-spy.co.uk/2012/05/install-rpi-gpio-python-library/>

TROUBLESHOOTING

White Bar Across Left Half of Screen

If a white bar appears across the left side of the screen and cannot click on the menu then follow these steps:

Control + Alt + F1

username: pi, password: ICP

`sudo rm -r /.config/lxpanel`

`startx`

HDMI Will Not Display

The Raspberry Pi must first be plugged in before the HDMI can be plugged in.

Code Will Not Execute Because Pi Cannot Discover I2C Bus Location

Chances are that the DAC chip has fried and needs to be replaced. If this is not the case, check to see which I2C location the Pi is trying to output the code by typing into the terminal "i2cdetect -y 1".

Pi Will Not Turn On

Check to see what color is being displayed by the LED. If only red is showing, the SD card is either inserted incorrectly or corrupted. It should have a red LED with a green LED that intermittently blinks.

Board Goes Off Intermittently

This is due to the power. Make sure that the Pi is plugged into a power supply that outputs 5V and 2.5A.

DIGITAL ANALOG CONVERTING SOFTWARE

The MCP4725 is a 12-bit DAC chip that converts the digital signal into an analog output to be properly displayed on the ProPAQ patient monitor. Determining what output voltage is given is. The voltage that is output is determined by the following equation:

$$V_{out} = \frac{(V_{cc} * Bit \ Value)}{4096}$$

The V_{cc} is the output voltage from the Raspberry Pi. The system has a V_{cc} of 5V. There are two sources from the Raspberry Pi, 3.3V and 5V. The 5V source output higher resolution waveform and that is why that source was chosen.

FREQUENCY CALCULATION

In order for the ICP waveform to have the same time synchronization as the EKG waveform, they must have the same frequency. This is calculated using the following equation:

$$Length \ of \ Waveform = e^{\frac{-(BPM+193.74)}{54.91}}$$

This equation can be seen in each of the waveform model files (ex. modelone.py).

COMPUTER APPLICATION

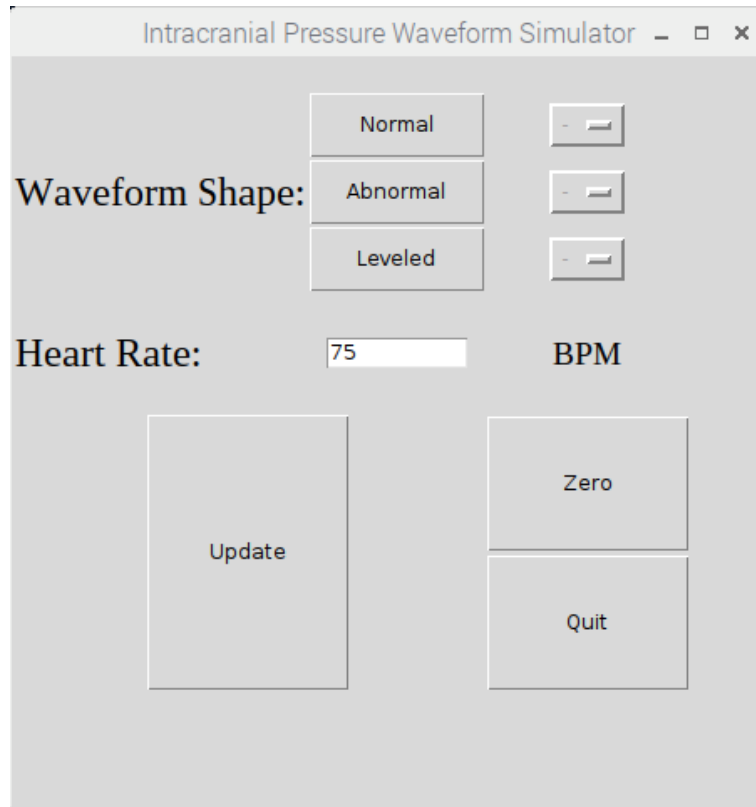


Figure 1: GUI

Listing 1: Gui Code

```
1 from guizero import Combo, App, Box, Text, Slider, TextBox, ButtonGroup, PushButton,
   CheckBox, MenuBar
2 import subprocess
3 import signal
4 import os
5 import sys
6
7 #Pulled in the necessary tools to create the GUI
8
9 class MyGUI:
10     def __init__(self, master):
11         self.master = master
12         #Creates the application window
13
14         self.space2 = Text(master, text=" ", grid=[0,0], width="fill")
15
16
17         self.Shapemessage1 = Text(master, size=18, text="Waveform Shape:", grid
=[0,2], align="left", font="Times")
18
19         self.Waveform1 = PushButton(master, text="Normal", grid=[1,1], width = "10",
command = self.show_choices1)
20         self.Waveform2 = PushButton(master, text="Abnormal", grid=[1,2], width="10",
command = self.show_choices2)
21         self.Waveform3 = PushButton(master, text="Leveled", grid=[1,3], width="10",
command = self.show_choices3)
22
23         #self.Waveform1.toggle()
24         #self.Waveform2.toggle()
25         #self.Waveform3.toggle()
26
27         self.combo1 = Combo(master, options=["-", "10", "11", "12", "13", "14", "15", "16", "
17"], grid=[2,1], selected = "-")
28         self.combo2 = Combo(master, options=["-", "21", "22", "23", "24", "25", "26"], grid
=[2,2])
29         self.combo3 = Combo(master, options=["-", "0"], grid=[2,3])
30
31         self.combo1.disable()
32         self.combo2.disable()
33         self.combo3.disable()
34         #Creates the three waveform buttons as well as the selection for ICP value
```

```
35     #When the waveform selected the user is then able to select an ICP value
36
37     self.space2 = Text(master, text=" ", grid=[0,6])
38     #Creates space between Waveform Shape and Heart Rate
39
40     self.HRmessage = Text(master, size=18, text="Heart Rate:", grid=[0,10], align
41     ="left", font="Times")
42     #Creates the HR subtitle
43
44     self.rate = TextBox(master, grid=[1,10])
45     self.HRunit = Text(master, size=15, text="BPM", grid=[2,10], font="Times")
46     #Creates the HR textbox to allow for typed input
47
48     self.space3 = Text(master, text=" ", grid=[0,11], font="Times")
49     #Creates space between Heart Rate and the Update Button
50
51     self.button = PushButton(master, command=self.update_value, width="12",
52     height="10", text="Update", grid=[0,13,2,2], align="bottom")
53     #Creates an Update Button that updates the inputed information and exports
54     the user-inputs to the next part of the program
55     self.when_clicked = self.update_value
56     #self.when_mouse_enters = self.update_value
57     #Creates the update button, no new code will be executed until update button
58     is pressed
59
60     self.button = PushButton(master, command=self.calibration, width="12", height
61     ="4", text="Zero", grid=[2,13], align="bottom")
62     #Creates zero button which is used to calibrate
63
64     self.button = PushButton(master, command=self.quit, width="12", height="4",
65     text="Quit", grid=[2,14], align="bottom")
66     #Creates the quit button which will exit out of app and stop the process from
67     running
68
69     self.pid = -1 #Initialize process id
70
71     #Launches the function for modelone , Normal ICP, in a seperate python process
72     and obtains the process ID
73     def modelone(self):
74         self.done()
```

```
69 #Launches the function for modeltwo , Leveled ICP, in a seperate python process
and obtains the process ID
70 def modeltwo(self):
71     self.done()
72
73 #Launches the function for modelthree, Abnormal ICP, in a seperate python process
and obtains the process ID
74 def modelthree(self):
75     self.done()
76
77 #Launches the calibration
78 def calibration(self):
79     self.done()
80     self.process = subprocess.Popen('python calibration.py', shell=True,
preexec_fn=os.setsid)
81     self.pid = self.process.pid
82
83 #Used to show the ICP value choices for the normal waveform
84 def show_choices1(self):
85     self.done()
86     if self.Waveform1.value == 0:
87         self.combo1.enable()
88         self.combo2.disable()
89         self.combo3.disable()
90     if self.Waveform1.value == 1:
91         self.combo1.disable()
92         self.combo1.value = 0
93
94 #Used to show the ICP value choices for the abnormal waveform
95 def show_choices2(self):
96     self.done()
97     if self.Waveform2.value == 0:
98         self.combo2.enable()
99         self.combo1.disable()
100         self.combo3.disable()
101     if self.Waveform2.value == 1:
102         self.combo2.disable()
103         self.combo2.value = 0
104
105 #Used to show the ICP value choices for the leveled waveform
106 def show_choices3(self):
107     self.done()
```

```
108         if self.Waveform3.value == 0:
109             self.combo3.enable()
110             self.combo1.disable()
111             self.combo2.disable()
112         if self.Waveform3.value == 1:
113             self.combo3.disable()
114             self.combo3.value = 0
115
116     #Command for updating the heart rate value
117     def update_value(self):
118         self.done()
119         if self.rate.value == "":
120             self.rate.value=65
121         if self.combo1.value == "10":
122             self.process = subprocess.Popen('python modelone.py {}'.format(self.rate.
value), shell=True, preexec_fn=os.setsid)
123             self.pid = self.process.pid
124         if self.combo1.value == "11":
125             self.process = subprocess.Popen('python modelone.py {}'.format(self.rate.
value), shell=True, preexec_fn=os.setsid)
126             self.pid = self.process.pid
127         if self.combo1.value == "12":
128             self.process = subprocess.Popen('python modelone.py {}'.format(self.rate.
value), shell=True, preexec_fn=os.setsid)
129             self.pid = self.process.pid
130         if self.combo1.value == "13":
131             self.process = subprocess.Popen('python modelone.py {}'.format(self.rate.
value), shell=True, preexec_fn=os.setsid)
132             self.pid = self.process.pid
133         if self.combo1.value == "14":
134             self.process = subprocess.Popen('python modelone.py {}'.format(self.rate.
value), shell=True, preexec_fn=os.setsid)
135             self.pid = self.process.pid
136         if self.combo1.value == "15":
137             self.process = subprocess.Popen('python modelone.py {}'.format(self.rate.
value), shell=True, preexec_fn=os.setsid)
138             self.pid = self.process.pid
139         if self.combo1.value == "16":
140             self.process = subprocess.Popen('python modelone.py {}'.format(self.rate.
value), shell=True, preexec_fn=os.setsid)
141             self.pid = self.process.pid
142         if self.combo1.value == "17":
```

```
143         self.process = subprocess.Popen('python modelone.py {}'.format(self.rate.  
value), shell=True, preexec_fn=os.setsid)  
144         self.pid = self.process.pid  
145  
146  
147  
148  
149  
150         if self.combo2.value == "21":  
151             self.process = subprocess.Popen('python modelthree.py {}'.format(self.  
rate.value), shell=True, preexec_fn=os.setsid)  
152             self.pid = self.process.pid  
153         if self.combo2.value == "22":  
154             self.process = subprocess.Popen('python modelthree.py {}'.format(self.  
rate.value), shell=True, preexec_fn=os.setsid)  
155             self.pid = self.process.pid  
156         if self.combo2.value == "23":  
157             self.process = subprocess.Popen('python modelthree.py {}'.format(self.  
rate.value), shell=True, preexec_fn=os.setsid)  
158             self.pid = self.process.pid  
159         if self.combo2.value == "24":  
160             self.process = subprocess.Popen('python modelthree.py {}'.format(self.  
rate.value), shell=True, preexec_fn=os.setsid)  
161             self.pid = self.process.pid  
162         if self.combo2.value == "25":  
163             self.process = subprocess.Popen('python modelthree.py {}'.format(self.  
rate.value), shell=True, preexec_fn=os.setsid)  
164             self.pid = self.process.pid  
165         if self.combo2.value == "26":  
166             self.process = subprocess.Popen('python modelthree.py {}'.format(self.  
rate.value), shell=True, preexec_fn=os.setsid)  
167             self.pid = self.process.pid  
168  
169  
170         if self.combo3.value == "0":  
171             self.process = subprocess.Popen('python modeltwo.py {}'.format(self.rate.  
value), shell=True, preexec_fn=os.setsid)  
172             self.pid = self.process.pid  
173  
174  
175         self.combo1.disable()  
176         self.combo2.disable()
```

```
177         self.combo3.disable()
178
179
180         self.combo1.value = "-"
181         self.combo2.value = "-"
182         self.combo3.value = "-"
183
184
185     #Kills the process that is currently running
186     def done(self):
187         if self.pid > 0:
188             try:
189                 os.killpg(os.getpgid(self.pid), signal.SIGTERM)
190                 self.pid = -1
191             except ProcessLookupError:
192                 pass
193     def quit(self):
194         self.done()
195         self.master.destroy()
196 root = App(title="Intracranial Pressure Waveform Simulator", width=450, height=450,
197           layout="grid")
198 #Pulls up the GUI
199 my_gui = MyGUI(root)
200 root.display()
201 #Creates then application window
```

WAVEFORM CODES


The waveform codes each rely on a different MATLAB file which contains a graph of the specified waveform. It then converts each point on the MATLAB plot into bit values that the DAC will recognize and be able to output the specified voltage.

The Calibration Code is used to zero the ProPAQ. This should be done before running the simulation just to initialize the system. This should be run before zeroing the probe on the ProPAQ. The process is as follows:

1. Start the GUI and press the zero button
 2. Zero the probe on the ProPAQ and change the scaling from 30 to -30
- execute the desired waveform, ICP value, and heart rate

Listing 2: Calibration

```
1 #Imports to Recognize the Dac on the Raspberry Pi
2 import Adafruit_MCP4725
3
4 dac = Adafruit_MCP4725.MCP4725()
5
6
7 zero = 300;
8 zero = int(zero)
9 dac.set_voltage(zero) #Set the output Voltage to Zero
```

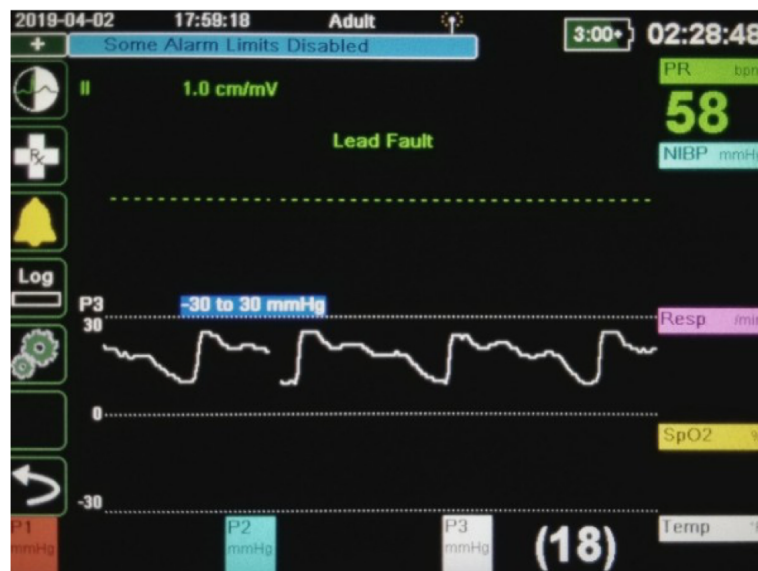**Figure 2: Normal Waveform**

Listing 3: Normal Waveform

```
1 from scipy import io
2 import numpy as np
3 import time
4 import sys
5 import math
6 import Adafruit_MCP4725
7
8 dac = Adafruit_MCP4725.MCP4725()
9 first_arg = sys.argv[1]
10
11 def modelone(intvalue = first_arg): #first_arg
12     intvalue = int(intvalue)
13     mat = io.loadmat('modelonescale.mat');
14     y_value = mat['interpoaltdatay'];
15     y_value = np.transpose(y_value);
16     new_value = np.zeros(y_value.shape)
17
18     for jj in range(len(y_value)):
19         if y_value[jj] == 28:
20             new_value[jj] = 208
21         elif y_value[jj] == 27:
22             new_value[jj] = 206
23         elif y_value[jj] == 26:
24             new_value[jj] = 204
25         elif y_value[jj] == 25:
26             new_value[jj] = 202
27         elif y_value[jj] == 24:
28             new_value[jj] = 200
29         elif y_value[jj] == 23:
30             new_value[jj] = 198
31         elif y_value[jj] == 22:
32             new_value[jj] = 196
33         elif y_value[jj] == 21:
34             new_value[jj] = 194
35         elif y_value[jj] == 20:
36             new_value[jj] = 192
37         elif y_value[jj] == 19:
38             new_value[jj] = 189
39         elif y_value[jj] == 18:
40             new_value[jj] = 187
41         elif y_value[jj] == 17:
```

```
42         new_value[jj] = 184
43     elif y_value[jj] == 16:
44         new_value[jj] = 182
45     elif y_value[jj] == 15:
46         new_value[jj] = 179
47     elif y_value[jj] == 14:
48         new_value[jj] = 177
49     elif y_value[jj] == 13:
50         new_value[jj] = 174
51     elif y_value[jj] == 12:
52         new_value[jj] = 172
53     elif y_value[jj] == 11:
54         new_value[jj] = 169
55     elif y_value[jj] == 10:
56         new_value[jj] = 167
57     elif y_value[jj] == 9:
58         new_value[jj] = 164
59     elif y_value[jj] == 8:
60         new_value[jj] = 162
61     elif y_value[jj] == 7:
62         new_value[jj] = 158
63     elif y_value[jj] == 6:
64         new_value[jj] = 156
65     elif y_value[jj] == 5:
66         new_value[jj] = 154
67     elif y_value[jj] == 4:
68         new_value[jj] = 152
69     elif y_value[jj] == 3:
70         new_value[jj] = 149
71     elif y_value[jj] == 2:
72         new_value[jj] = 147
73     elif y_value[jj] == 1:
74         new_value[jj] = 144
75
76     var = 1;
77
78     X = math.exp(-(intvalue+193.74)/54.91)
79     new_value = new_value.flatten()
80     while var == 1:
81         for val in new_value:
82             val = int(val)
83             dac.set_voltage(val)
```

```
84         time.sleep(X)
85     if __name__ == "__main__":
86         modelone()
```

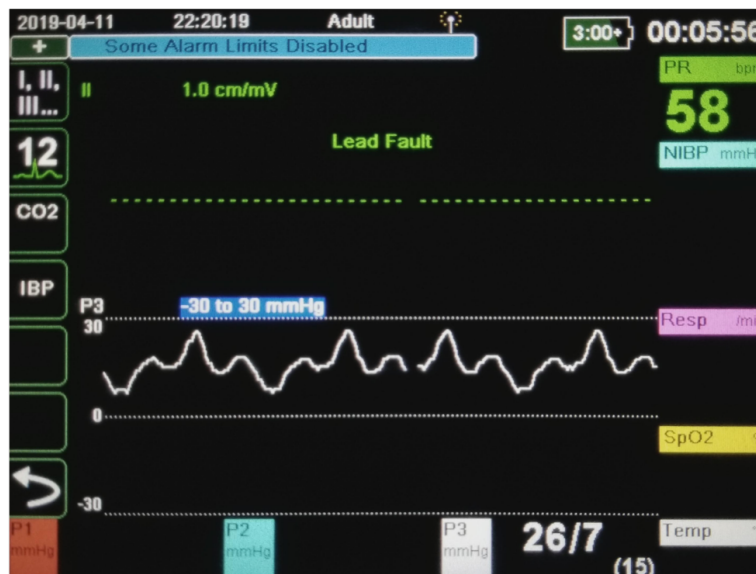


Figure 3: Abnormal Waveform

Listing 4: Abnormal Waveform

```
1 from scipy import io
2 import numpy as np
3 import time
4 import sys
5 import math
6 import Adafruit_MCP4725
7
8 dac = Adafruit_MCP4725.MCP4725()
9 first_arg = sys.argv[1]
10
11 def modelthree(intvalue = first_arg):
12     intvalue = int(intvalue)
13     mat = io.loadmat('ICPL.mat');
14     y_value = mat['ICPL'];
15     y_value = np.transpose(y_value);
16     new_value = np.zeros(y_value.shape)
17
18     for jj in range(len(y_value)):
19         if y_value[jj] == 31:
20             new_value[jj] = 219
21         elif y_value[jj] == 30:
22             new_value[jj] = 217
23         elif y_value[jj] == 29:
24             new_value[jj] = 214
25         elif y_value[jj] == 28:
26             new_value[jj] = 212
27         elif y_value[jj] == 27:
28             new_value[jj] = 202
29         elif y_value[jj] == 26:
30             new_value[jj] = 199
31         elif y_value[jj] == 25:
32             new_value[jj] = 196
33         elif y_value[jj] == 24:
34             new_value[jj] = 184
35         elif y_value[jj] == 23:
36             new_value[jj] = 182
37         elif y_value[jj] == 22:
38             new_value[jj] = 179
39         elif y_value[jj] == 21:
40             new_value[jj] = 177
41         elif y_value[jj] == 20:
```

```
42         new_value[jj] = 167
43     elif y_value[jj] == 19:
44         new_value[jj] = 164
45     elif y_value[jj] == 18:
46         new_value[jj] = 162
47     elif y_value[jj] == 17:
48         new_value[jj] = 160
49     elif y_value[jj] == 16:
50         new_value[jj] = 150
51     elif y_value[jj] == 15:
52         new_value[jj] = 148
53     elif y_value[jj] == 14:
54         new_value[jj] = 146
55     elif y_value[jj] == 13:
56         new_value[jj] = 144
57     elif y_value[jj] == 12:
58         new_value[jj] = 142
59     elif y_value[jj] == 11:
60         new_value[jj] = 139
61     elif y_value[jj] == 10:
62         new_value[jj] = 138
63     elif y_value[jj] == 9:
64         new_value[jj] = 136
65     elif y_value[jj] == 8:
66         new_value[jj] = 134
67     elif y_value[jj] == 7:
68         new_value[jj] = 132
69     elif y_value[jj] == 6:
70         new_value[jj] = 130
71     elif y_value[jj] == 5:
72         new_value[jj] = 128
73     elif y_value[jj] == 4:
74         new_value[jj] = 126
75     elif y_value[jj] == 3:
76         new_value[jj] = 120
77     elif y_value[jj] == 2:
78         new_value[jj] = 116
79     elif y_value[jj] == 1:
80         new_value[jj] = 114
81
82 new_value = new_value.flatten()
83 var = 3;
```

```

84 X = math.exp(-(intvalue+193.74)/54.91)
85 while var == 3:
86     for val in new_value:
87         val = int(val)
88         dac.set_voltage(val)
89         time.sleep(X)
90
91 if __name__ == "__main__":
92     modelthree()

```

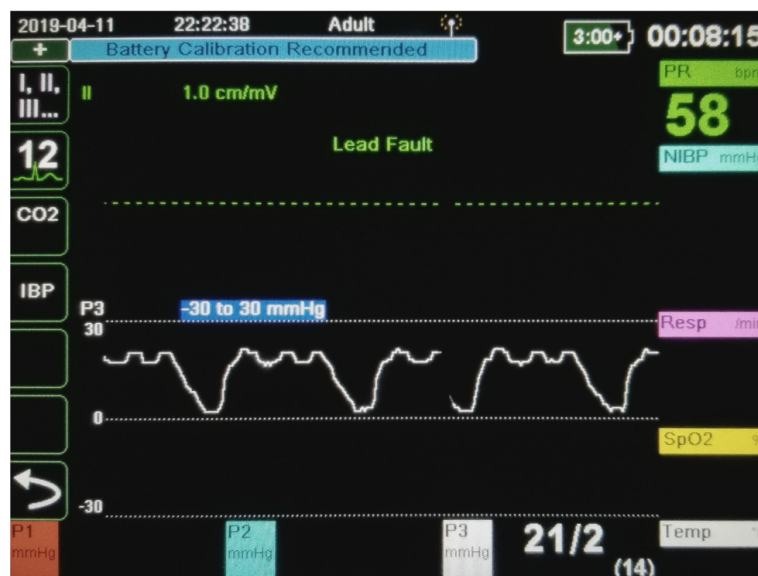


Figure 4: Leveled Waveform

Listing 5: Leveled Waveform

```
1 from scipy import io
2 import numpy as np
3 import time
4 import sys
5 import math
6 import Adafruit_MCP4725
7
8 dac = Adafruit_MCP4725.MCP4725()
9 first_arg = sys.argv[1]
10
11 def modeltwo(intvalue = first_arg):
12     intvalue = int(intvalue)
13     mat = io.loadmat('modeltwoscale.mat');
14     y_value = mat['modeltwo'];
15     y_value = np.transpose(y_value);
16     new_value = np.zeros(y_value.shape)
17
18     for jj in range(len(y_value)):
19         if y_value[jj] == 31:
20             new_value[jj] = 219
21         elif y_value[jj] == 30:
22             new_value[jj] = 217
23         elif y_value[jj] == 29:
24             new_value[jj] = 214
25         elif y_value[jj] == 28:
26             new_value[jj] = 212
27         elif y_value[jj] == 27:
28             new_value[jj] = 202
29         elif y_value[jj] == 26:
30             new_value[jj] = 199
31         elif y_value[jj] == 25:
32             new_value[jj] = 196
33         elif y_value[jj] == 24:
34             new_value[jj] = 184
35         elif y_value[jj] == 23:
36             new_value[jj] = 182
37         elif y_value[jj] == 22:
38             new_value[jj] = 179
39         elif y_value[jj] == 21:
40             new_value[jj] = 177
41         elif y_value[jj] == 20:
```

```
42         new_value[jj] = 167
43     elif y_value[jj] == 19:
44         new_value[jj] = 164
45     elif y_value[jj] == 18:
46         new_value[jj] = 162
47     elif y_value[jj] == 17:
48         new_value[jj] = 160
49     elif y_value[jj] == 16:
50         new_value[jj] = 150
51     elif y_value[jj] == 15:
52         new_value[jj] = 148
53     elif y_value[jj] == 14:
54         new_value[jj] = 146
55     elif y_value[jj] == 13:
56         new_value[jj] = 144
57     elif y_value[jj] == 12:
58         new_value[jj] = 142
59     elif y_value[jj] == 11:
60         new_value[jj] = 139
61     elif y_value[jj] == 10:
62         new_value[jj] = 138
63     elif y_value[jj] == 9:
64         new_value[jj] = 136
65     elif y_value[jj] == 8:
66         new_value[jj] = 134
67     elif y_value[jj] == 7:
68         new_value[jj] = 132
69     elif y_value[jj] == 6:
70         new_value[jj] = 130
71     elif y_value[jj] == 5:
72         new_value[jj] = 128
73     elif y_value[jj] == 4:
74         new_value[jj] = 126
75     elif y_value[jj] == 3:
76         new_value[jj] = 120
77     elif y_value[jj] == 2:
78         new_value[jj] = 116
79     elif y_value[jj] == 1:
80         new_value[jj] = 114
81
82 new_value = new_value.flatten()
83 var = 2;
```



```
84     X = math.exp(-(intvalue+193.74)/54.91)
85     while var == 2:
86         for val in new_value:
87             val = int(val)
88             dac.set_voltage(val)
89             time.sleep(X)
90 if __name__ == "__main__":
91     modeltwo()
```